This Page Is Inserted by IFW Operations and is not a part of the Official Record

BEST AVAILABLE IMAGES

Defective images within this document are accurate representations of the original documents submitted by the applicant.

Defects in the images may include (but are not limited to):

- BLACK BORDERS
- TEXT CUT OFF AT TOP, BOTTOM OR SIDES
- FADED TEXT
- ILLEGIBLE TEXT
- SKEWED/SLANTED IMAGES
- COLORED PHOTOS
- BLACK OR VERY BLACK AND WHITE DARK PHOTOS
- GRAY SCALE DOCUMENTS

IMAGES ARE BEST AVAILABLE COPY.

As rescanning documents will not correct images, please do not report the images to the Image Problem Mailbox.

THIS PAGE BLANK (USPTQ)



WORLD INTELLECTUAL PROPERTY ORGANIZATION International Bureau



D3(A)

INTERNATIONAL APPLICATION PUBLISHED UNDER THE PATENT COOPERATION TREATY (PCT)

(51) International Patent Classification 5:

G01N 1/28

A1

(11) International Publication Number: WO 91/18273

(43) International Publication Date: 28 November 1991 (28.11.91)

(21) International Application Number:

PCT/US91/03389

(22) International Filing Date:

15 May 1991 (15.05.91)

(30) Priority data:

524,339

16 May 1990 (16.05.90)

US

(71) Applicant: SCIENTIFIC IMAGING INSTRUMENTS, INC. [US/US]; 30 Lindeman Drive, Trumbull, CT 06611 (US).

(72) Inventor: AUSTIN, I., Mark; 31 Pierce Street, New Rochelle, NY 10801 (US).


(74) Agents: GREASON, Edward, W. et al.; Kenyon & Kenyon, One Broadway, New York, NY 10004 (US).

(81) Designated States: AT, AT (European patent), AU, BB, BE (European patent), BF (OAPI patent), BG, BJ (OAPI patent), BR, CA, CF (OAPI patent), CG (OAPI patent), CH, CH (European patent), CI (OAPI patent), CM (OAPI patent), DE, DE (European patent), DK, DK (European patent), ES, ES (European patent), FI, FR (European patent), GA (OAPI patent), GB, GB (European patent), GR (European patent), HU, IT (European patent), JP, KP, KR, LK, LU, LU (European patent), MC, MG, ML (OAPI patent), MR (OAPI patent), MW, NL, NL (European patent), NO, PL, RO, SD, SE, SE (European patent), SN (OAPI patent), TG (OAPI patent), TG (OAPI patent).

Published

With international search report.

(54) Title: METHOD AND APPARATUS FOR PREPARING A LIQUID SPECIMEN

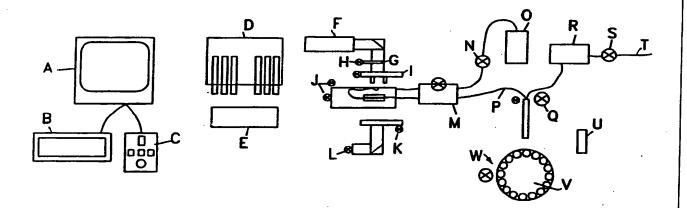

(57) Abstract

A method and apparatus are provided for accurately preparing centrifuged fluid samples for microscopic analysis. In the method a centrifuged fluid specimen contained in a sample tube of known volume and dimensions located in motor (W) driven sample tube cassette carousel (V) has the level of fluid specimen in the sample tube detected therein with a sensing duel pipette tube (D6). This level is used in calculating the total volume of fluid specimen in the sample tube with the computer (D). Thereafter, an amount of fluid specimen sufficient to leave a predetermined volume percentage of fluid specimen in the sample tube is decanted. Next the sample tube is agitated to resuspend the sediment that had been concentrated in the centrifuged fluid specimen.

FOR THE PURPOSES OF INFORMATION ONLY

Codes used to identify States party to the PCT on the front pages of pamphlets publishing international applications under the PCT.

AT	Austria	50			
AU	Australia	ES	Spain	MG	Madagascar
BB	Barbados	. Pi	Finland	ML	Mali
BE	Belgium	FR	France	MN	
BF	•	GA	Gabon	MR	Mongolia
BG	Burkina Faso	GB	United Kingdom		Mauritania
	Bulgaria	GN	Guinea	MW	Malawi
BJ	Benin	GR	Greece	NL	Netherlands
BR	Brazil	HU	Hungary	, NO	Norway
CA	Canada	n		PL	Poland
CF	Central African Republic	JP	ltaly	RO	Romania
CG	Congo	_	Japan	SD	Sudan
,CH	Switzerland	KP	Democratic People's Republic	SE	Sweden
CI	Côte d'Ivoire		of Korea	SN	Senegal
CM	Cameroon	KR	Republic of Korea	SU	Soviet Union
CS	Czechosłovakia	Li	Licchtenstein	TD	
DE		LK	Sri Lanka	TG	Chad
DK	Germany	LU	Luxembourg		Togo
UK	Denmark	MC	Моласо	us	United States of America

METHOD AND APPARATUS FOR PREPARING A LIQUID SPECIMEN

BACKGROUND OF THE INVENTION

5 The present invention relates generally to the field of automated instrumentation, and more specifically to the field of microscopic analysis.

One known system for microscopically analyzing fluids is disclosed in United States Letters Patent No. 4,804,267 to Greenfield, assigned to the assignee of the present invention, the disclosure of which is incorporated herein by reference. Although the Greenfield system provided effective solutions to many of the problems confronting the art it itself possesses several disadvantages and drawbacks.

Specifically, the system of the '267 patent utilizes a single pump and flushes the specimen through the system to waste. Experience has shown that this may allow for the introduction of bubbles into the flow cell which can be seen as artifacts under high magnification. In addition, the system of the '267 patent provides only a limited number of user functions or features.

There is a long felt need, which has gone unsatisfied prior to the making of the present invention, for an automated system for microscopically analyzing fluids,

20

which provides rapid, accurate, multiple sample viewing capability, and the ability to control the viewing and analysis of the specimen including saving it, if necessary. Also there exists an unsatisfied, long felt need for a method and means for accurately preparing a fluid specimen for microscopic analysis.

It is accordingly a general object of the present invention to overcome the aforementioned limitations and drawbacks associated with known systems and to fulfill the needs mentioned above by providing a system for microscopically analyzing a specimen having all of the desirable attributes noted above.

It is a particular object of the present invention to provide a method for accurately preparing fluid specimens for microscopic analysis.

It is another object of the present invention to provide an automated method for preparing fluid specimens for microscopic analysis.

Another object of the present invention is to provide a method for standardizing the preparation of fluid specimens for microscopic analysis.

A further object of the present invention is to provide a means of mixing (resuspending) the fluid sediment.

The foregoing and other objects and advantages which
will be apparent in the following detailed description
of the preferred embodiment, or in the practice of the
invention, are achieved by the invention disclosed
herein, which generally may be characterized as a method
for preparing for microscopic analysis a centrifuged
fluid specimen contained in a sample tube of known
volume and dimensions comprising the steps of:

WO 91/18273 PCT/US91/03389

-3-

detecting the level of fluid specimen in the sample tube; calculating the total volume of fluid specimen in the sample tube; decanting an amount of fluid specimen to leave a predetermined volume percentage of fluid 5 specimen remaining in the sample tube; and agitating the predetermined percentage of fluid specimen remaining in the sample tube to resuspend the sediment that had been concentrated in the centrifuged fluid specimen, and apparatus for preparing for microscopic analysis a 10 centrifuged fluid specimen contained in a sample tube of known volume and dimensions comprising: means for detecting the level of fluid specimen in the sample tube; means for calculating the total volume of fluid specimen in the sample tube; means for decanting an 15 amount of fluid specimen to leave a predetermined volume percentage of fluid specimen in a sample tube; and means for agitating the predetermined percentage of fluid specimen remaining in the sample tube to resuspend the sediment that had been concentrated in the centrifuged 20 fluid specimen.

BRIEF DESCRIPTION OF THE DRAWINGS

Serving to illustrate exemplary embodiments of the invention are the drawings, of which:

Figure I is a block diagram of the robotic microscope of the present invention;

Figure 2 is a diagram illustrating the internal components of the three subsystems of the robotic microscope of the present invention;

Figure 3 is an exploded diagram of a multi-channel flow 30 cell in accordance with the present invention;

Figure 4 is block diagram of a multi-channel flow cell with pumping/sampling system in accordance with the present invention;

Figure 5 is a block diagram of a multi-channel pipette
and pumping/sampling system in accordance with the present invention; and

Figure 6 is a block diagram of apparatus for preparing a specimen in accordance with the present invention.

DETAILED DESCRIPTION OF THE PREFERRED EMBODIMENT

- As used herein, and in accordance with the present invention, a robotic microscope is a computerized optical imaging instrument which automates various aspects of microscopic laboratory tests and is designed to assist the trained laboratory technician in
- performing routine microscopic analysis, such as, for example, urinalysis. Among the benefits over conventional approaches are improved accuracy and standardization; improved ability to visualize difficult specimens; reduced specimen handling; reduced
- 20 disposables; and higher productivity for the laboratory.

Referring to Figure 1, a diagram of the three subsystems of the robotic microscope of the present invention is illustrated. As shown therein, the robotic microscope consists of a Control Station, an Imaging Station and a Specimen Station.

Referring now to Figure 2, the internal components of the three subsystems of the robotic microscope are illustrated. As shown therein, the Control Station includes a monitor for viewing the microscopic specimen, 30 a keyboard and other computer input device (such as a mouse or trackball) as necessary to interface with the

computer controls of the microscope. The Control Station is where the laboratory technologist controls the microscope, performs the analysis and records the results.

- The screen on the monitor displays a highly magnified image of the specimen as well as relevant data about the specimen and alphanumeric and graphical symbols which, when used in conjunction with the proper input device, allow a technician to control all the parameters
- required for the microscopic examination. In particular, while the magnified specimen is being viewed, the technician can also control the following functions: magnification; focus; location; scanning; sampling; optical enhancement; lighting; and sample data entry.

The Imaging Station includes a high magnification optical and video system, a computer, input/output interface, power supply, filters, flow cell, apertures, necessary motors and control electronics. The Imaging Station is where a magnified electronic image of the specimen to be analyzed is obtained.

The Specimen Station includes apparatus to deliver the specimen to the viewing area of the Imaging Station and means for removing the specimen from the Imaging Station after it has been analyzed. In the instance of solid or dried specimen on a slide, the Specimen Station consists of a mechanical device to load and unload slides. In the instance where the specimen is wet or in liquid state, the Specimen Station includes a pumping system to pump specimen into and flush specimen from the viewing area of the Imaging Station and an indexing system to index sequential samples.

The Specimen Station also includes a carousel which holds the specimen cassette tray. Each cassette holds up to 20 specimens in sample tubes and is coded and numbered to assure accurate transcription of patient information. Directly behind the cassette is the washbowl, and above both is the sampling autopipette. To the left of the carousel and autopipette is the flush reservoir, which is accessed via the top cover door.

The majority of all operations performed on the robotic microscope are done with the Trackball and the corresponding cursor on the Main Control Screen.

The operation of the robotic microscope illustrated in Figures 1 and 2 is as follows. After turning on the power switch and waiting for instrument to complete its automatic priming and self-testing functions, the technician loads the sample tubes into the Sample Tube Cassette (V). The technician then uses the trackball (C) to activate the computer controls on the monitor screen (A) and instructs the instrument to prepare the next sample.

The instrument then indexes the Carousel Motor (W) until Sample Tube Cassette (V) has turned a sample tube into position under the Two Axis motorized Pipette (Q). The pipette (Q) then decants the specimen and pumps it via tubing (P), valves (M) through the flow cell (J) via the sampling pump (N).

The technician then performs the analysis of the specimen, viewing it on screen (A) via computer (D), Camera (F), Optical Column (G), Lens (I), flow cell (J), Apertures (K), and Illuminator (L).

After completing the analysis, the technician enters the report via keyboard (B), and the system is automatically

flushed via sample pump (N), Flush Reservoir (0), out to washbowl (U). The waste pump (S) then takes the specimen and flush solution from the washbowl (U) through the pipette (Q) out through the waste port (T). The instrument is now ready for the next sample.

All of these functions are controlled by the technician through the computer (D) and powered by the power supply (E).

In accordance with the present invention, multiple
sample viewing capability for liquid samples is provided
by means of a multi-channel flow cell. Although the
following discussion of the multi-channel flow cell is
in terms of a dual channel flow cell, it is clear that
the number of channels can be increased accordingly.

- 15 The general operation of the dual channel flow cell is as follows. The specimen is first pumped into the alpha channel of the flow cell. While the operator is examining the specimen under optical magnification on the monitor screen the pump flushes and then loads the beta channel of the flow cell. When the operator has completed due examination of the specimen contained in the alpha channel, the Imaging Station moves the beta channel of the flow cell into view under the optical system. This provides the operator with extremely rapid 25 access to sequential prepared samples because the time required to prepare the second sample is coincident with the operator's time to view the previous sample, thus there is no waiting time for the operator while the sample is prepared.
- Referring to Figures 3 and 4, an exploded diagram of a dual channel flow cell, and a dual channel flow cell with pumping/sampling system, respectively, are illustrated. As shown therein, the dual channel flow

€...

cell requires an input line and an output line for both alpha and beta channels and a valve system on each end to switch the pump between channels. The valves also guarantee that the specimen will be held rigidly in place during microscopic examination.

The dual channel flow cell consists of two sample channels or chambers the depth of which are determined by the working distance and top cover thickness given the optical characteristics of the objectives. The thickness of the top cover is also determined by the objective while the thickness of the bottom is determined by the working distance of the condenser lens.

The dual channel flow cell includes transparent upper and lower retaining members, generally flat in form, one 15 of which has a plurality of pairs of fluid flow passages formed therein. It also includes a central member, generally flat in form, having a plurality of display chambers defined therein. Each of the display chambers is in fluid flow registration with one of the pairs of fluid flow passages and each of the display chambers is out of fluid flow registration with all of the other of the plurality of display chambers. The upper, lower and central members are secured to one another to form an integral structure, which is carried in a body having a 25 flat central well and a viewing aperture formed therein. At least a portion of each of the plurality of display chambers underlies the viewing aperture.

The operation of the multi-channel flow cell is as

follows. (1) The sample is taken up in pipette D4

through the first channel of the two channel pipette via

the action of sample pump E4. (2) After the sample has

passed through valve C4 and into the alpha channel of

the Flow cell A4, the action of pump E4 ceases and then

the valve C4 switches to the other channel (beta). The pipette is then inserted into the next sample and steps 1 and 2 are repeated with the exception that the sample is pumped into the beta channel and valve C4 subsequently switches to channel alpha. (4) To remove the sample, the pipette is inserted into the washbowl (not shown). (5) The sample pump is then operated in reverse mode and pumps flush solution out of flush solution reservoir G4 through the alpha channel of flow 10 cell A4, through valve C4 and through first channel of pipette D4, thereby pushing sample into washbowl and filling alpha channel of the flow cell with flush solution. (6) The waste pump F4 is then activated and pulls specimen and flush solution out of the washbowl up 15 through the second channel of pipette D4 and out through waste port H4. (7) Steps 4 through 6 are repeated to clean out the beta channel of the flow cell with the exception that the valve C4 is first set to the beta (8) The same procedures as listed above can also be used for more than two channel flow cells, requiring only more channels in the flow cell itself and greater switching capacity in the valve.

In accordance with the present invention, a means of sampling liquid specimens such that each sample channel of a multi-channel pipette and its associated tubing acts as a reservoir of specimen is provided.

Referring to Figure 5, a multi-channel pipette and pumping/sampling system in accordance with the present invention is illustrated. As shown therein, the specimen is sampled via a pipette which descends into the sample tube. The pipette consists of three rigid tubes or channels, two for sample and an auxiliary one for waste. Each sample channel in the pipette corresponds to one of the channels in the flow cell.

35 Flexible tubing of appropriate diameter connects the

rigid tubes of the pipette to the valves, and then from the valves to the flow cell.

The volume of specimen contained in the rigid and flexible tubing functions as a reservoir, holding additional sample. This allows more sample to be viewed in addition to the sample that is in the viewing channel of the flow cell. Further, this reservoir of specimen in the tubing allows the sample to be saved in its entirety as opposed to being lost in the flushing 10 These design advantages are juxtaposed to a system that has only a single tube for taking up the Such a system would require a valve to switch between the two channels of the flow cell. advantages of speed from the two channels of the flow cell were to be maintained, then the balance of the specimen would either be left in the sample tube (requiring a separate operation to extract it if additional specimens were to be viewed) or the specimen would be lost.

- The multi-channel pipette system, described in the system shown in Figure 5 has all of the same functional capabilities as the multi-channel flow cell system of Figure 4, but provides additional capabilities as well. In particular, by dedicating one channel of the pipette directly to one channel of the flow cell, the tubing can function as a reservoir of additional sample thereby providing additional system features such as sample advance; sample saving; sample staining all of which require the sample reservoir to be functional.
- The operation of the multi-channel pipette is as follows. The system functions the same as the system described in Figure 4 above with the exception that a separate sample channel in the pipette and its associated tubing correspond to one of the sample

channels in the multi-channel flow cell. In addition, it has the following operations. (1) To advance the sample into the flow cell, the valve C5 is switched to the appropriate channel of the flow cell being viewed. (2) The sample pump is then run in the direction to pull the sample towards the pump by an amount equal to the volume of the channel of the flow cell. (3) The valve is then switched back to the other channel. (4) Sample saving is similar to the flushing operation described in the system of Figure 4, with the exception that the pipette is placed back in the original sample tube instead of the washbowl, and the waste pump is never activated. (5) Sample staining is the same as sample saving with the exception that after the sample has been saved, the technician then adds stain to the sample and 15 then it is resampled to the flow cell in the procedure described in Figure 4 above.

1

In accordance with the present invention, the system also provides for the capability to rapidly and 20 accurately prepare a suspension of sample from a centrifuged fluid specimen, such as a urine specimen. After the specimen is centrifuged, the biological sediment is concentrated in the bottom of the sample The dual or multi-channel pipette then descends 25 into the sample tube, and, utilizing a sample channel and the auxiliary (waste) channel in the pipette as sensor probes, detects the level of the fluid in the sample tube and determines the amount of fluid in the sample tube. The two channels consist of rigid metal 30 tubing electrically isolated from each other within the pipette and are used with known circuitry to detect the changes in conductivity between the air and the fluid in the sample tube. The computer then calculates the total volume of fluid in the sample tube and determines the amount of fluid to be decanted in order to leave a 35 predetermined percentage (for example 10%) in the sample

tube. The pipette then descends to this predetermined level and decants the fluid through the auxiliary channel in the pipette along the way. The pipette then descends further into the sample tube and the auxiliary channel in the pipette and the waste pump then cycle vigorously. This back and forth cycling of the pump agitates the button of biological sediment in the bottom of the sample tube creating a concentrated suspension. This system allows for an exact suspension to be prepared independently of the initial volume presented. In practice, the current manual approaches are extremely lax in the precision with which the concentration is prepared.

Referring to Figure 6, a method and apparatus for

15 accurately preparing liquid samples that have
centrifuged sediment in them in accordance with the
present invention is illustrated. Although a multichannel flow cell is shown therein, it is noted that the
system does not require the use of a multi-channel flow

20 cell.

Its operation is as follows: (1) The pipette D6
descends into the sample tube (not shown) and detects
the level of sample liquid. (2) The system calculates
the volume and determines the amount to be decanted to

25 create a proportional specimen (i.e. a 10% suspension of
sediment and original sample supernatant). The system
decants this proportion via Valve F6 and pump G6. (3)
The system then cycles via valve F6 and Pump G6 (that is
pumps rapidly and forcefully backwards and forwards)

30 with the pipette immersed in the sample. This breaks up
centrifuged constituents and causes them to create a
suspension of said constituents that is highly and
proportionally concentrated in comparison to the
original liquid volume. (4) The specimen is then

WO 91/18273 PCT/US91/03389

-13-

sampled as in a procedure described for either system in Figure 4 or 5.

Although the present invention has been described in terms of its presently preferred embodiment, certain modifications thereof based on the descriptions and teachings herein may be apparent to those skilled in the art. For example, the embodiment disclosed above deals with a system for microscopically analyzing fluids. An adaptation of the invention to other forms of specimens such as solid samples should be apparent to those skilled in the art.

Accordingly, the scope of the present invention is defined by the following claims.

processors that the contract of the first

WHAT IS CLAIMED IS:

1. A method for preparing for microscopic analysis a centrifuged fluid specimen contained in a sample tube of known volume and dimensions comprising the steps of:

detecting the level of fluid specimen in the sample tube;

calculating the total volume of fluid specimen in the sample tube;

decanting an amount of fluid specimen to leave a predetermined volume percentage of fluid specimen in the sample tube; and

agitating the predetermined percentage of fluid specimen remaining in the sample tube to resuspend the sediment that had been concentrated in the centrifuged fluid specimen.

- A method as recited in Claim 1 wherein the steps thereof are automated.
- 3. Apparatus for preparing for microscopic analysis a centrifuged fluid specimen contained in a sample tube of known volume and dimensions comprising:

means for detecting the level of fluid specimen in the sample tube;

means for calculating the total volume of fluid specimen in the sample tube;


means for decanting an amount of fluid specimen to leave a predetermined percentage volume of fluid

specimen in the sample tube; and

means for agitating the predetermined percentage of the fluid specimen remaining in the sample tube to resuspend the sediment that had been concentrated in the centrifuged fluid specimen.

- 4. Apparatus as recited in Claim 3 wherein said detecting means, said calculating means, said decanting means and said agitating means are controlled by a computer.
- 5. Apparatus as recited in Claim 4 wherein said detecting means includes pipette means configured to detect changes in conductivity between the air and the fluid.
- 6. Apparatus as recited in Claim 5 wherein said calculating means includes computer means.
- 7. Apparatus as recited in Claim 6 wherein said decanting means include pipette means having at least an auxiliary channel, waste pump means and waste valve means, said auxiliary channel of said pipette means being connected to said waste pump means and waste valve means by fluid flow connection means.
- 8. Apparatus as recited in Claim 7 wherein said agitating means include said pipette means having at least an auxiliary channel, said waste pump means and said waste valve means, said auxiliary channel of said pipette means being connected to said waste pump means and said waste valve means by fluid flow connection means.

1/6

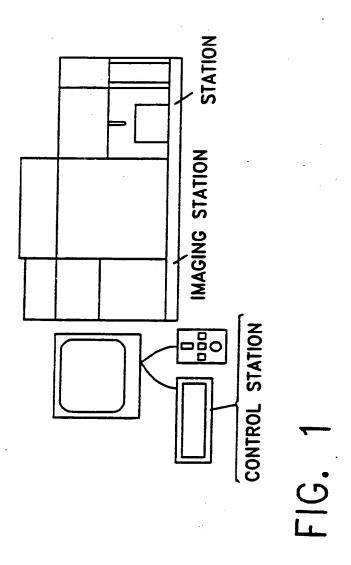

SUBSTITUTE SHEET

2/6



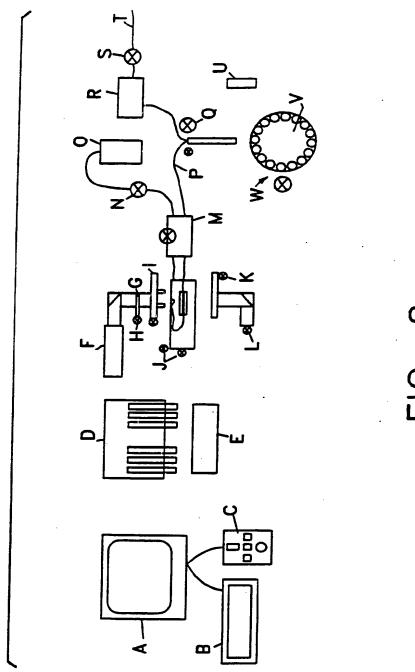


FIG. 2


3/6

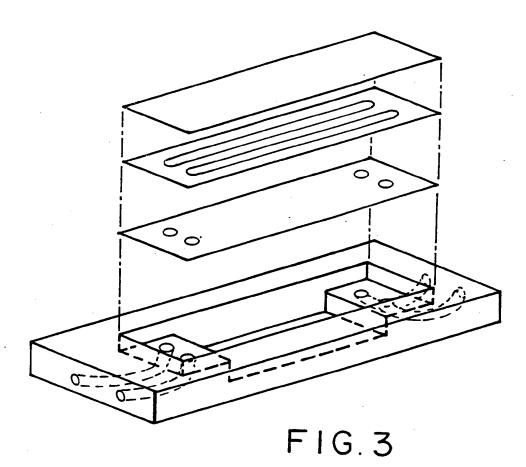

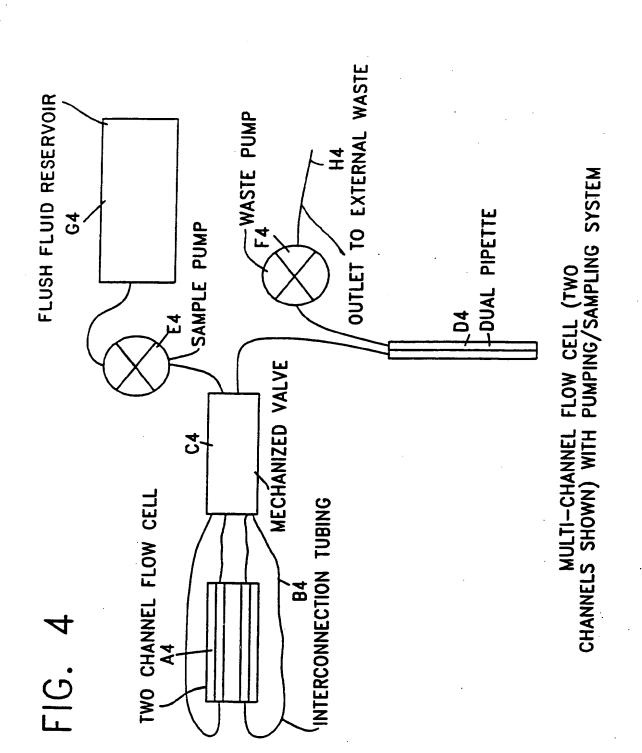

5/6

FLUSH FLUID RESERVOIR OUTLET TO EXTERNAL WASTE MULTI-CHANNEL PIPETTE F5 WASTE PUMP H5 **G**2 SAMPLE PUMP MULTI – CHANNEL PIPETTE AND PUMPING SAMPLING SYSTEM Ŕ MECHANIZED VALVE S TWO CHANNEL FLOW CELL NTERCONNECTION TUBING **B**2 FIG. 5 **A**5

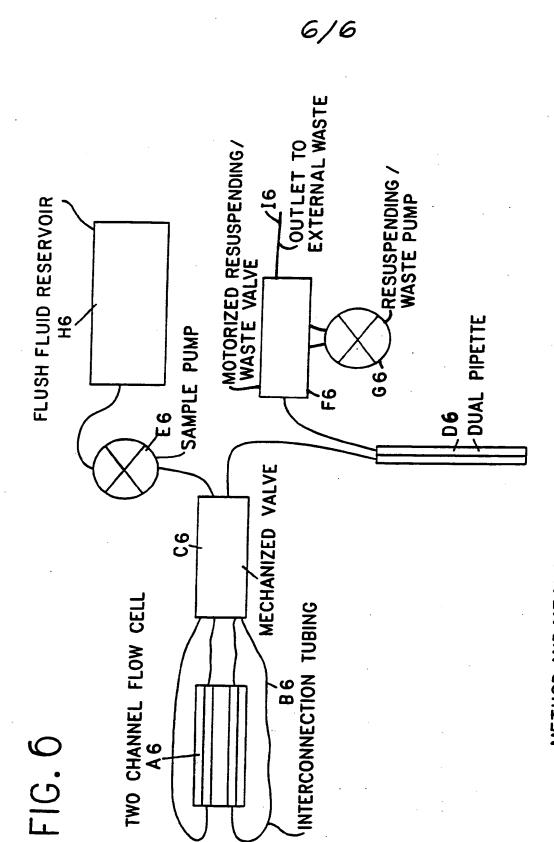

METHOD AND MEANS FOR PREPARING SPECIMEN

INTERNATIONAL SEARCH REPORT

International Application No. PCT/US91/03389 I. CLASSIFICATION OF SUBJECT MATTER (if several classification symbols apply, indicate all) 6 US CL: 73/863,803.21,863.32,863.02;356/36 IPC(5): GO1N 1/28II. FIELDS SEARCHED Minimum Documentation Searched 7 Classification System Classification Symbols 73/863,863.21,863.32,863.01,863.02,864.11,864.15,864.21, U.S. 864.22; 356/36, 366/140; 436/45,174,177; 422/64 Documentation Searched other than Minimum Documentation to the Extent that such Documents are Included in the Fields Searched III. DOCUMENTS CONSIDERED TO BE RELEVANT 9 Category • Citation of Document, 11 with indication, where appropriate, of the relevant passages 12 Relevant to Claim No. 13 US, A, 4,435,293 (GRAHAM JR., ET AL.) 06 MARCH 1984 Y 1-6 Note col. 1, lines 14-65 and col. 9, line 56-col. 11, line 24. US, A, 4,436,631 (GRAHAM JR., ET AL.) 13 MARCH 1984, Note col. 1 lines 17-56; col. 8, lines 19-38 and col. Y 1-6 12, lines 1-36. Y US, A, 4,486,315 (TEIPEL) 04 DECEMBER 1984, Note col. 6, line 58-col. 7, line 5. 1-6 US, A, 4,939,925 (SAKUMA ET AL.) 10 JULY 1990 Note abstract, Fig. 1, and Fig. 2. Y,P 1-6 Special categories of cited documents: 10 later document published after the international filing date or priority date and not in conflict with the application but cited to understand the principle or theory underlying the invention document defining the general state of the art which is not considered to be of particular relevance earlier document but published on or after the international filing date "X" document of particular relevance; the claimed invention cannot be considered novel or cannot be considered to involve an inventive step document which may throw doubts on priority claim(s) or which is cited to establish the publication date of another citation or other special reason (as specified) document of particular relevance; the claimed invention cannot be considered to involve an inventive step when the document is combined with one or more other such documents, such combination being obvious to a person skilled document referring to an oral disclosure, use, exhibition or ments, su in the art. document published prior to the international filing date but later than the priority date claimed "4" document member of the same patent family IV. CERTIFICATION Date of the Actual Completion of the International Search 1 6 AUG 1991 18 MAY 1991 International Searching Authority Signature of Authorized Officer

TOM NOLAND

ISA/US

m PCT/ISA/210 (second sheet) (Rev.11-67)

PCT/US	91/03389
SHEET)	

ategory *	Citation of Document, with indication, where appropriate, of the relevant passages	
Y	the relevant passages	Relevant to Claim No
	Patent Abstracts of Japan Grp. p +78, vol.13, no. 572, 18 December 1989, (abstract of 01-240859,) See entire document.	1-6
Y !	JP, A, 01-240859, (Sakuma et al.) 26 September 1989 See fig. 1 and fig. 2.	1-6
Y	IIS, A, 4,873,633 (MEZEI ET AL.) 10 OCTOBER 1989, Note abstract, Fig. 1, Fig. 2, Fig. 5A, Fig. 5B, col. 14, line 51-col. 15, line 18, and col. 26, lines 47-53.	1-6
A	US, A, 4,693,972 (MANSOUR ET AL.) 15 SEPTEMBER 1987.	1,3
A	US, A, 3,754,444 (URE ET AL.) 28 AUGUST 1973	2,5,6
A	US, A, 4,487,836 (TAKAYANAGI ET AL.) 11 DECEMBER 1984.	5
A	Patent Abstracts of Japan Grp. P. 699, vol. 12, no. 151, 11 May 1988 (Abstract of 62-269037).	5
4	JP, A, 62-269037 (Nakahara) 21 November 1987	5
1	US, A, 4,829,837 (Telfer) 16 May 1989	5
-		
	·	
		•
		1
		•

THIS PAGE BLANK (USPTO)

XP-002175195

AN - 1984-053958 [09]

AP - SE19820003349 19820601

CPY - VANR-I

DC - D16 J04 P81 S03 S05

FS - CPI;GMPI;EPI

IC - C12M1/36; G02B21/28

IN - VAN R

MC - D05-H D05-H02 J04-B01

- S03-E04B1 S03-E14H S05-C

PA - (VANR-I) VAN R

PN - SE8203349 A 19840116 DW198409 017pp

PR - SE19820003349 19820601

XA - C1984-023019

XIC - C12M-001/36; G02B-021/28

XP - N1984-040593

- AB SE8203349 Incubator comprises a closed container with a space for the specimen. The specimen can be illuminated for examination. The container can be used with different optical arrangements. The microscope can be used either inside or outside of the container.
 - The container can be kept at a constant temperature by a heating element. A gas pipe introduces air with an increased level of carbondioxide. There are two water chambers to provide a constant moisture level in the contaner. There is a door at each end. The temp. is monitored by a sensor and current supplied to the heater. The microscope objective is placed adjacent to the bottom of the container. The microscope is mounted on top of the container. The container is clamped to the table by clamps.
 - Incubator is versatile and can be used, e.g. in medicine. (Provisional Basic previously advised in Week 8405)(0/9)
- IW CLOSE CONTAINER MICROSCOPE ACT INCUBATE OBJECT STUDY ALLOW MICROSCOPE PART INSERT HELD NO LOSS SEAL
- IKW CLOSE CONTAINER MICROSCOPE ACT INCUBATE OBJECT STUDY ALLOW MICROSCOPE PART INSERT HELD NO LOSS SEAL

INW - VAN R

NC - 001

OPD - 1982-06-01

ORD - 1984-01-16

PAW - (VANR-I) VAN R

 TI - Closed container used with microscope - acts as incubator for object studied, and allows microscope parts to be inserted or held outside it with no loss of seal THIS PAGE BLANK (USPTO)